

EXCELIGHT®

Take your business to the
next level in advanced
treatments...

585 & 650 IPL



Lynton

EXCELIGHT® IPL

Made in Britain by the UK's No.1 IPL manufacturer & leading NHS supplier, the Lynton EXCELIGHT® offers businesses the opportunity to introduce the UK's most in-demand advanced aesthetic treatments in one compact and easy to use device.

Offering permanent hair reduction and a wide range of skin rejuvenation treatments, the Lynton EXCELIGHT® IPL is widely regarded by many Dermatologists as the most clinically effective IPL in the marketplace. Delivering industry-leading power and speed alongside easy to use, built-in treatment parameters, the Lynton EXCELIGHT® has been designed to ensure all practitioners can operate competently and carry out safe and effective IPL treatments.

Two highly effective IPL handpieces offering multiple treatment options

Designed to fit straight into any business looking to introduce advanced treatments, promising rapid Return on Investment

Performing over 10 different treatments, including permanent hair reduction, vascular lesions, pigmentation & skin rejuvenation



“

I have never looked back after purchasing the IPL from Lynton. It has helped my business expand from a very small salon with three treatment rooms, to new premises with nine treatment rooms. Lynton's training and service is excellent. We have achieved fantastic results and have many happy clients. I would recommend Lynton to anyone wanting to start with a hair removal or skin rejuvenation device.

Becky Parkin



INTERFACE PRO^a




Intuitive pre-set parameters result in a quick learning curve for all operators. The touch-screen treatment controls (above) make the EXCELIGHT® simple and easy to use, which helps all practitioners feel confident in delivering advanced IPL treatments.

Manual override is available for all experienced practitioners wishing to tailor the output of the system for an individual client.

Lynton Intense Pulsed Light (IPL)

The 585 and 650 handpieces are supplied as standard and are the perfect combination for hair removal and skin rejuvenation treatments.

Handpieces are easy to change over, which, combined with the simple touchscreen interface, allows you to quickly select the treatment and protocols you require.

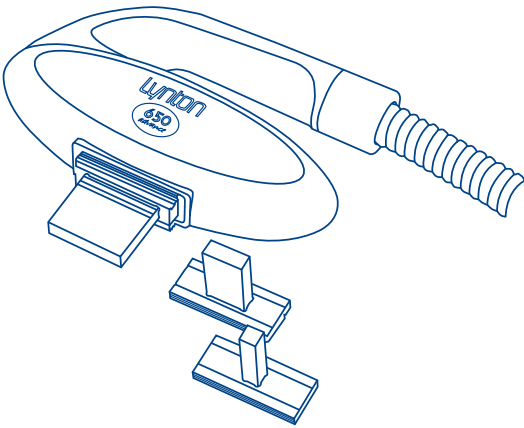

Applications with the EXCELIGHT® IPL:

- Hair Removal
- Thread Veins
- Rosacea
- Age Spots
- Skin Rejuvenation
- Skin Tightening
- Pigmentation
- Port Wine Stains

+ more



ILG^a

ILG™: Interchangeable Light Guide technology provides you with speed, precision and enhanced epidermal protection from bulk-heating. ILG™ allows the practitioner to tailor the spot size for every client due to its easy to use, snap-in design. Unrivalled spatially-selective targeting reduces discomfort in treatments, allowing more effective fluences to be safely used. ILG™ also provides effortless positioning in contoured areas, such as a small red vein on the client's nose, or a large spot size for hair removal on a client's back.



Pigmentation

Courtesy of The Lynton Clinic



Before

After

Hair Removal

Courtesy of Loretta Redding: Erazor Clinic



Before

After

Rejuvenation

Courtesy of The Court House Clinic



Before

After

Red Vein Removal

Courtesy of Kerry Belba | Laser Skin Solutions



Before

After

To learn more, or to arrange a demonstration, call Lynton on 01477 536 977 or visit lynton.co.uk



Call Lynton today on 01477 536 977
or visit lynton.co.uk

System Specifications

EXCELIGHT®

IPL (Intense Pulsed Light) HP	585	650
Wavelength (nm)	550 - 1100 nm	625 - 1100 nm
Pulse Train Energy (J)	10 - 250 J	10 - 250 J
Individual Pulse Durations (IPL)	<1ms – 10 ms	<1ms – 10 ms
Pulse Sequencing	2 – 6 pulses	2 - 6 pulses
Pulse Delay (ms)	10 – 60 ms	10 – 60 ms
Repetition Rate (Hz)	Up to 1 Hz	Up to 1 Hz
Spot Size (cm)	5.0 x 1.0, 2.5 x 1.0, 1.0 x 1.0 (standard) & 1.0cm Circular (optional)	
Fluence (J/cm ²)	Up to 40 J/cm ²	Up to 36 J/cm ²

Physical Parameters

Dimensions (cm)	37 (w) x 65 (d) x 49 (h)
Handpiece	Interchangeable Quartz Blocks
Control Panel	Full colour LCD and Touchscreen
Weight (kg)	Approx 40
Manufacturer	Lynton, UK

The Medical Devices Directive (Medical CE-marking) and Notified Body approvals for the EXCELIGHT Intense Pulsed Light System only relate to the medical indications detailed in the accompanying Operator Manual.

Lynton Lasers Ltd

Lynton House, Manor Lane,
Holmes Chapel, Cheshire, CW4 8AF

Tel: +44 (0)1477 536 977

info@lynton.co.uk

www.lynton.co.uk

October 2025



0123