




MOTUS AY **Recognize Your Star Player**



Focusing on Effectiveness and Comfort
with Moveo

Motus stands out in **Aesthetics** and **Dermatology**

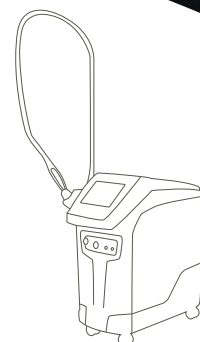

DEKA
Innate Ability



MOTUS AY Quick Wins:

Hair Removal | Vascular Lesions
Bening Pigmented Lesions | Onychomycosis
Non ablative Skin Rejuvenation
Pseudofolliculitis Barbae (PFB) | Leg Veins

Key Principles

Motus AY uses a combination of the most common wavelengths: **Alex 755 nm** – the most effective in hair removal and **Nd:YAG 1064 nm** in vascular treatments.

Why Choosing Motus:

Double Wavelength (Alex + Nd:YAG)

Wide range of applications: not only hair removal and pigmented lesions but also vascular treatments, to meet the growing demand of patients.

Moveo Technology

New frontier in hair removal and pigmented lesions: guarantee of safety, comfort, speed and results.

Incomparable Quality to Price Ratio

A succesfull combination of effectiveness and cost and low maintenance.



THE SCIENCE BEHIND



100% of Energy Absorption with Moveo Technology

The new handpiece technology integrates a sapphire tip that comes in contact with the skin, reducing the reflected energy loss. The new hair removal method consists in repeated passes on treated areas that cause gradual heating of the vital parts of hairs.

Effectiveness + Safety

The absence of energy dispersion means lower power for the same result.

Speed + Comfort

Multipass action allows to treat big areas in few minutes.

New plug-in
connection system



Moveo PL: The New Moveo Handpiece for Pigmented Lesions


Max results are obtained by exploiting the ideal wavelength in pigmentary (Alex 755 nm) in a **low fluence-multipass** delivery process with integrated cooling, giving the greatest comfort and safety to patients. MOVEO PL is safe, extremely effective and allows to treat any skin type in very few sessions, **with minimal side effects or downtime.**

NEW!
—
MOVEO PL



Choose up more than 10 Standard Handpieces

2.5-20 mm spots available with automatic spot recognition system and ergonomic connector for air-cooling.

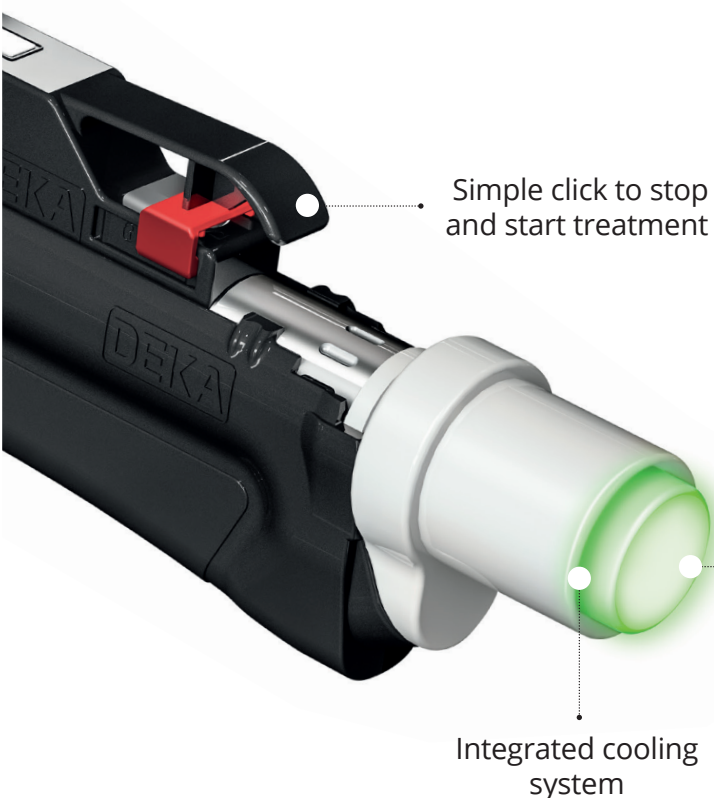

Simple click to stop
and start treatment

20 mm
sapphire tip

Integrated cooling
system

THE SCIENCE BEHIND



Revolutionary Technology Innovation

Motus AY makes hair removal with Alex effective, painless and simple for both **light phototypes** and **the darker ones**, as well as with thin and light hair.

The Nd:YAG source broadens the field of action by offering the possibility of dealing with **vascular lesions** such as capillaries, or telangiectasias.

Unique in-motion Alex and Nd:YAG lasers

Smoke free procedure

Double energy delivered to skin

Treatment for all phototypes



effectiveness



speed



safety




simplicity




Clinical Results



Alexandrite laser with Moveo Technology.
Courtesy of **M. Troiano**, M.D.
and **Prof. P. Bonan**, M.D. – Florence, Italy

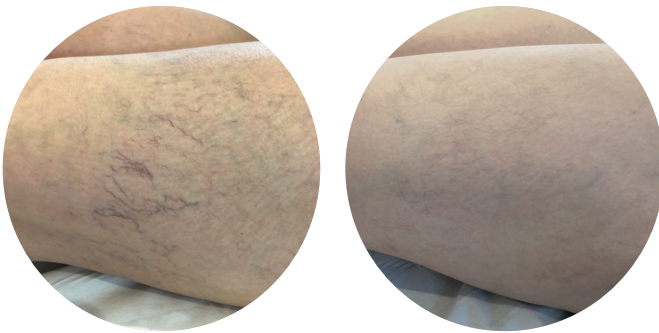

Nd:YAG laser source.
Courtesy of **Setsuko Beauty
Clinic Ltd**,
Bridgwater - United Kingdom

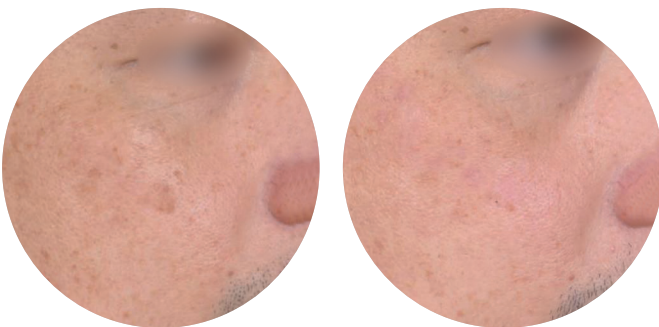

Moveo PL handpiece.
Courtesy of **Li Ying Lin**,
M,D. - Elite Clinic,
Taipei-Taiwan



Knowing Our Audiences and Market

Thanks to the growing beauty consciousness among consumers, worldwide, the global hair removal market is witnessing a high-paced increase in its size and valuation.

GLOBAL HAIR REMOVAL*:



1st
REQUESTED TREATMENT
WITHIN AESTHETIC INDUSTRY



+17%
COMPOUND ANNUAL
GROWTH RATE



OVER
3.9 billion \$
EXPECTED UNTIL 2026

*Sources: Market Study Report "Global Laser Hair Removal Market Size Status And Forecast 2019-2025"
Acumen Research and Consulting "Laser Hair Removal Market Size, Share, Trends, Scope, Growth Opportunity and Forecast 2019-2026"

The Floor to Practitioners

“ In the last several years we’ve had lasers that are safe for dark skin, but the limit has been efficacy. But the Motus AY was designed to reduce the laser’s energy loss as it moves over the skin. The concentrated, doubled energy transfer to skin means the Motus works better than other lasers on darker skin tones.”

Dr. Paul Jarrod Frank, Dermatologist,
New York – USA

“ With Motus AY DEKA has changed the whole concept of alexandrite laser. The effectiveness of this laser source is scientifically and clinically recognized for photoepilation. Thanks to the innovative MOVEO Technology, laser hair removal with alexandrite laser is now not only EFFECTIVE but also SAFE and PLEASANT, suitable for darker skin types too”

Prof. Paolo Bonan M.D.
Villa Donatello Clinic, Florence – Italy

“ I had always avoided undergoing treatments laser hair removal for the problems that are usually related to dark skin. I decided to try Moveo and stayed positively affected by both the result and the speed of treatment: it is much faster than waxing!”

Fiona May - World Champion Long Jump

THE PERFECT MATCH

MOTUS AY Strength

Operator and patient total comfort benefits:

- Fastest and only **in-motion Alex** and **Nd:YAG**
- **Fluididity** and lightness of the contact handpiece
- **Easy** to learn and use
- Effective on **all phototypes**
- **No energy loss**
- **Fastest** treatment: treats 100 cm² in approx 25 sec.
- **No burning** or skin discolouration
- **No pain**

 **Only for DEKA Users:**
a great support in medical practice
DEKA *Club*

Moveo Tecnología

2.5-20 mm spot sizes hanpieces

More than **200 protocols developed**
for Moveo and standard handpieces

Alex 755 nm
laser source

Nd:YAG 1064 nm
laser source

16 A electrical
supply



Technical Specifications

MOTUS AY		
Laser Source	Alexandrite 755 nm	Nd:YAG 1064 nm
Energy per Pulse	31 J (max)	63 J (max)
Laser Power	65 W (max)	95 W (max)
Repetition Rate	From 1 to 10 Hz	From 0.3 to 10 Hz
Spot Sizes (optional)	2.5 mm, 5 mm, 7 mm, 10 mm, 12 mm, 14 mm, 15 mm, 16 mm, 18 mm, 20 mm; Moveo HR, Moveo PL	
Aiming Beam	Green Diode Laser @532 nm < 4 mW	
Emission Control	Fingerswitch or Footswitch	
Handpiece Cooling	Integrated Skin Cooler	
Dimensions and Weight	90 (H) cm x 51 (W) cm x 83 (D) cm 90 Kg	
Electrical Requirements	230 Vac – single phase – 50/60 Hz – 16 A	

CAUTION - Visible and invisible laser radiation. Avoid eye or skin exposure to direct or scattered radiation. Class 4 laser product.
This brochure is not intended for the market of USA.



Dealer stamp



DEKA M.E.L.A. s.r.l.

Via Baldanzese, 17 - 50041 Calenzano (FI) - Italy
Tel. +39 055 8874942 - Fax +39 055 8832884

DEKA Innate Ability

A spin-off of the El.En. Group, DEKA is a world-class leader in the design and manufacture of lasers and light sources for applications in the medical field. DEKA markets its devices in more than 80 countries throughout an extensive network of international distributors as well as direct offices in Italy, France, Japan and USA. Excellence is the hallmark of DEKA's experience and recognition garnered in the sphere of R&D in over thirty years of activity. Quality, innovation and technological excellence place DEKA and its products in a unique and distinguished position in the global arena. DEKA manufactures laser devices in compliance with the specifications of Directive 93/42/EEC and its quality assurance system is in accordance with the ISO 9001 and ISO 13485 standards.