




promax lipo

The ProMax Lipo utilises a powerful combination of three market-leading technologies for instant inch loss, dramatic cellulite reduction and non-surgical facial skin lifting & tightening.

Pain-free and highly profitable, the ProMax Lipo delivers fast treatment times with no consumables. This all-in-one platform helps give your clients maximum results on a variety of concerns, creating a superb return on investment opportunity.

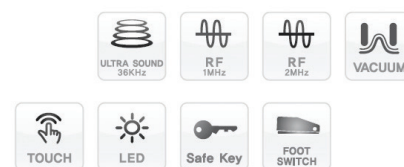

Call Lynton today on 01477 536 977
or visit lynton.co.uk

System Specifications

Ultrasonic Frequency	36 kHz
Radio Frequency	1MHz (Multi-Polar), 2MHz (Bipolar)
Vacuum Range	650 mmHg
Screen Panel	7 inch, Digital Touch Screen
Input Voltage	220-240V / 50-60Hz
Weight	52 kg
Mains supply	220 – 240 VAC, 13 A

Lynton Lasers Ltd

Lynton House, Manor Lane,
Holmes Chapel, Cheshire, CW4 8AF
Tel: +44 (0)1477 536 977
Fax: +44 (0)1477 536 978
info@lynton.co.uk

MAR 19 - SK

Aesthetics Awards
EQUIPMENT SUPPLIER OF THE YEAR



promax lipo

Maximise your body treatment results

Bipolar & Multipolar Radio Frequency, Vacuum Massage, Ultrasonic Cavitation & LED Technology



instant **inch loss**
dramatic **cellulite reduction**
body & facial **skin tightening**

Bipolar Radio Frequency & Vacuum for instant tightening & lifting on the face and around the eye and jawline areas.

Multi-Polar RF & Vacuum for firmer, tighter skin on the body, helping reduce the appearance of cellulite.

Large touchscreen interface with easy to use controls - designed to be safely delivered by all practitioners, even those with little or no experience.

Ultrasonic Cavitation tackles stubborn areas of fat that prove unresponsive to exercise, such as abdomen, upper arms and outer thighs.



The Science...



Ultrasonic Cavitation

Body

Sound waves penetrate deep into the fat layer causing vibrations and the formation of microbubbles within the interstitial fluid. As these bubbles collapse, the resulting shock wave leads to a temporary change in the fat cell membrane, allowing lipids to escape the cell and be removed by the body. Perfect for stubborn areas of fat that prove unresponsive to exercise, such as the abdomen, upper arms and outer thighs.

Multi-Polar Radio Frequency & Vacuum

Body


Uses controlled waves of energy to create heating in the subcutaneous and deep dermal layers. This thermal effect improves circulation to increase fat cell metabolism and reduce the appearance of cellulite. The thermal energy also results in instant collagen contraction and creates long term improvement to the collagen fibres for firmer, tighter skin. Vacuum massage technology provides powerful lymphatic drainage to process the lipids released during cavitation and to remove toxins to further improve the appearance of cellulite.

Bi-Polar Radio Frequency & Vacuum

Face

Creates thermal energy within facial dermal tissues for instant tightening and lifting. Perfect for around the eye, jawline and neck, not only does the ProMax Lipo instantly tighten, it also provides long term collagen stimulation for a younger, fresher appearance. The introduction of vacuum at this stage also improves circulation and removes toxins for a brighter complexion.

Main Treatment Areas

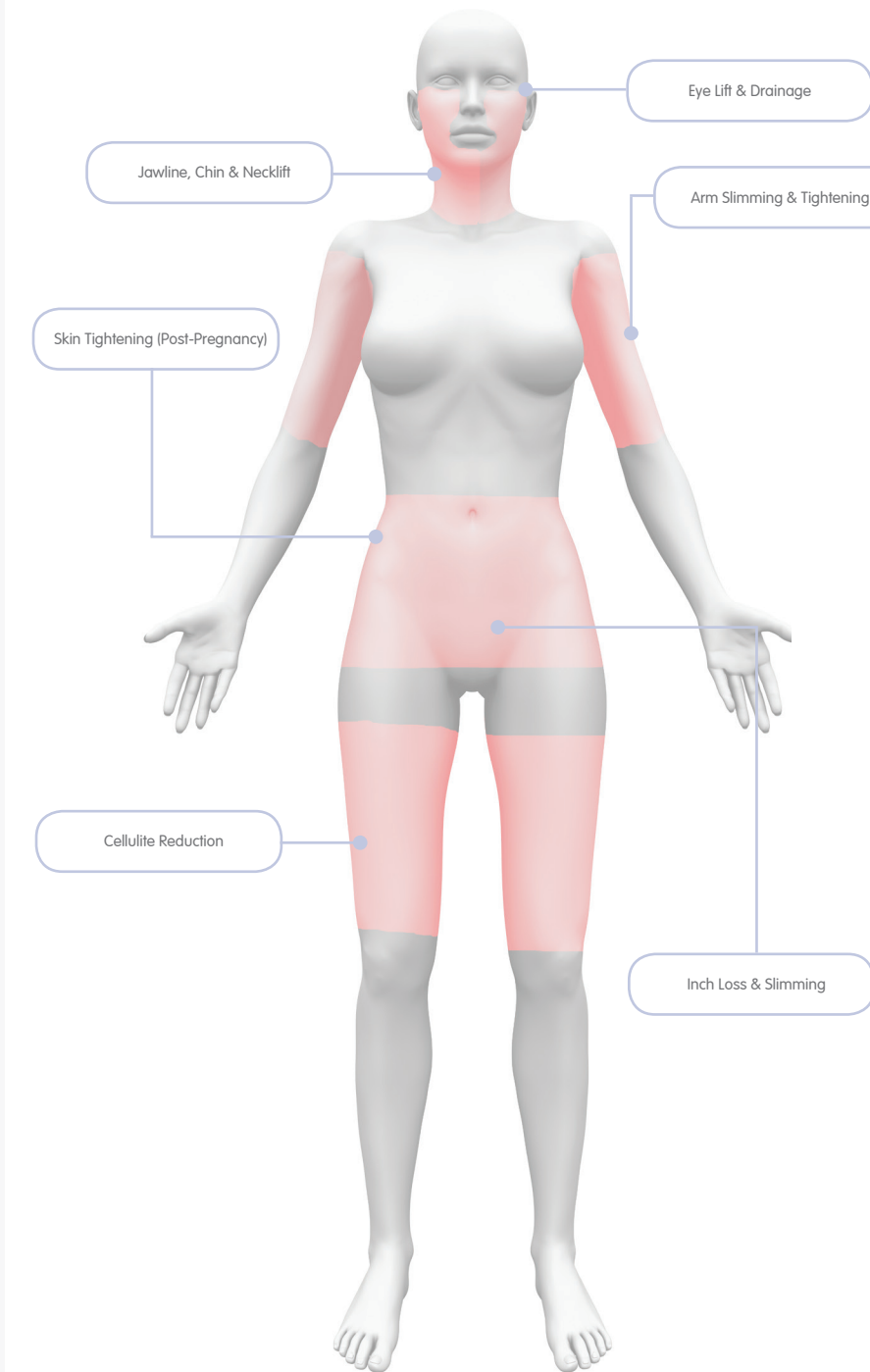


To learn more or arrange a demonstration, call Lynton on 01477 536 977 or visit lynton.co.uk/promax-lipo

Results

Body Shaping

5x Sessions | Courtesy of Laser Clinic North Wales

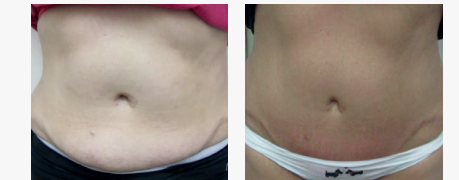

Body Shaping

5x Sessions | Courtesy of The Lynton Clinic



Body Shaping

5x Sessions | Courtesy of Luciderm



Skin Tightening (Body)

5x Sessions | Courtesy of Laser Clinic North Wales



Skin Tightening (Body)

6x Sessions | Courtesy of Laser Clinic North Wales



Skin Tightening (Face/Neck)

2x Sessions | Courtesy of The Lynton Clinic

