

DYNAMIX®

# Aesthetic Mixology

Plug & Play Your Way



Lynton

---

# DYNAMIX®

---

Dynamix is more than a device; it is a clinic in a box, tailor-made just for you to transform your business journey and fuel your clinic's growth.

This medically CE-marked device seamlessly blends IPL, RF, tattoo removal, hair removal, and advanced lasers for unparalleled results in a way that is unique to your clinic.

Two highly effective ReBright & ReMove handpieces, offering hair removal & **various skin rejuvenation treatment options**

ReMove (Long-pulsed Nd:YAG) available for those looking to treat **ALL skin types including darker skins**

**ReModel** radiofrequency causes micro-vibrations within the skin tissue to stimulate fibroblast cells which then produce new collagen and elastin, creating stronger, firmer dermal tissue.

Award-winning ResurFACE® handpiece upgrade, available for **fractional ablative treatments** with minimal downtime



Dynamix empowers you to become a hub of comprehensive aesthetic solutions, helping you to stand out in the competitive marketplace.

---

# Why Dynamix?

---

Dynamix is the key to unlocking a new era of success for your business:

- **British Craftsmanship:** Dynamix represents Lynton's 30 years of market-leading manufacturing, ensuring excellence in every single component.
- **All-in-One Efficiency:** Streamline your clinic operations with Dynamix's comprehensive suite of technologies - ReBright, ReMove, ReModel & ResurFACE® - all in a single device.
- **Excellent Return On Investment:** Compared to multiple devices, Dynamix offers an advanced solution, enabling you to maintain affordability without compromising on quality and results.



“

I knew I wanted to expand my aesthetic business and the multi-functional platform of the Dynamix by Lynton has enabled me to do exactly that! I started by offering facial treatments including fractional resurfacing, IPL for pigmentation, acne, and rosacea, as well as laser carbon peels to reduce pore size, all of which compliment my aesthetic business and regular patients.

As a nurse practitioner it was important to have a medical-grade device provided by a reputable company. The post sales and training support has been above and beyond my expectations from when I started this journey.

I am continually amazed by the compact size of the Dynamix considering how many hand pieces/treatments it offers. I would certainly recommend it as a safe and reliable option which provides outstanding results. ”

**Michelle Worthington**  
**MWAH - Michelle Worthington Aesthetic Health**  
**Sutton Coldfield, Birmingham**

To learn more or arrange a demonstration, call Lynton on 01477 536 977 or visit [lynton.co.uk](http://lynton.co.uk)



---

# Your Vision, Your Clinic

---

Your aesthetic clinic, your rules – Dynamix redefines versatility. It's not just a device, it's **your** device.

Dynamix is an all-in-one solution that is custom built for you to maximise your clinic space. No more juggling multiple devices or thinking which one to get; Dynamix sets the stage for your growth, and you can have it all, your way.

---

## ILG™

---

ILG™: Interchangeable Light Guide technology provides you with speed, precision and enhanced epidermal protection from bulk-heating.

ILG™ allows the practitioner to tailor the spot size for every client due to its easy-to-use, snap-in design. Unrivalled spatially-selective targeting reduces discomfort in treatments, allowing more effective fluences to be safely used.

ILG™ also provides effortless positioning in contoured areas, such as a small red vein on the client's nose, or a large spot size for hair removal on the client's back.



## HANDPIECE OPTIONS

### 1. ReBright



ReBright 585 IPL is Lynton's medical-grade IPL for delivering excellent and consistent clinical treatments, including the market-leading Illumifacial.

Standard applications with the ReBright 585 and ReMove 650:

- Acne
- Red Veins
- Pigmentation
- Skin Rejuvenation
- Sun Damage
- + more



### 2. ReMove

Your choice of both laser and IPL technologies for hair removal and skin rejuvenation on all skin types. Remove LP (Long-pulsed Nd:YAG) can also be used for skin laxity and vascular concerns.



### 3. ReMove Q

ReMove Passive Q Switch laser handpieces available in 532nm and 1064nm for fading multicolour tattoos and treating pigmentation.



### 4. ReModel®

ReModel radiofrequency causes micro-vibrations within the skin tissue to stimulate fibroblast cells which then produce new collagen and elastin, creating stronger, firmer dermal tissue. This leads to a subtle lifting and tightening of the skin that can greatly improve the loss of firmness.



### 5. ResurFACE®

Add this award-winning fractional Erbium YAG handpiece at 2940nm to achieve dramatic results on your clients. Use fractional ablative operating mode for a wide variety of applications, such as stretch marks, scarring, and fine lines & wrinkles.


---

# Dynamix Results

---

## Before & Afters

ReBright (585nm)  
IPL Skin Rejuvenation

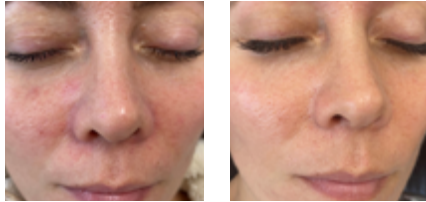

Before After  
Courtesy of ByNatashaTurner

ReMove (650nm)  
IPL Hair Removal




Before After  
Courtesy of Setsuko Beauty

ReMove (650nm)  
IPL Hair Removal




Before After  
Courtesy of Loretta Redding - Erazor Clinic

Remove 1064 Long Pulsed Nd YAG  
Leg Veins

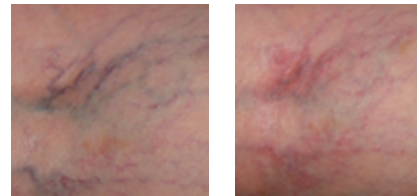

Before After  
Courtesy of Dr Lalloo, The Cosmetic & Medical Centre, Dublin

ResurFACE  
(Erbium YAG 2940)

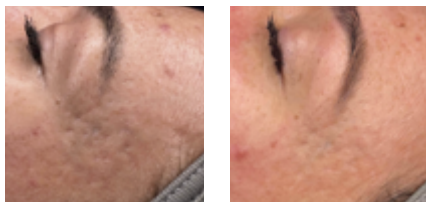

Before After  
Courtesy of The Skin Studio by Blush & Blow London

Remove Q 532 and 1064  
Tattoo Removal



Before After  
Courtesy of The Lynton Clinic

ReModel  
Radiofrequency



Before After  
Courtesy of The Lynton Clinic



# illumiFacial<sup>®</sup>

The illumiFacial<sup>®</sup> flips the script on the photofacial of the past and creates a new normal of medically-proven IPL to produce results with little/no downtime that will delight clients time after time. This completely non-invasive treatment uses a unique combination of specialist Lynton skincare which has been expertly formulated and tested for use alongside Lynton's medical-grade devices.

Designed to transform the skin in four steps the illumiFacial<sup>®</sup> offers an instant, long-lasting improvement for a variety of skin concerns:

- 1** Facial Cleanse  
The first step involves a deep cleanse with the Lynton Facial Cleanser, which removes all traces of dirt and impurities whilst protecting the lipid barrier.
- 2** Tri-Fruit Acid Peel  
Lynton's unique Tri-Fruit Acid Peel is applied to the skin for anywhere between 3 - 10 minutes. This peel resurfaces skin texture, decongests pores, and prepares the skin for the IPL treatment.
- 3** Lynton IPL  
Lynton's medical-grade IPL technology stimulates collagen and targets any discoloration on the skin. In most cases the results are instant, resulting in dramatically clearer, firmer skin.
- 4** Hydrate & Protect  
The IPL treatment is followed up with Lynton Light Soothe and Light Protect, specifically formulated to hydrate and protect the skin.





Call Lynton today on 01477 536 977  
or visit [lynton.co.uk](http://lynton.co.uk)

---

#### System Specifications

#### ReBright & ReMove

---

Wavelengths	ReBright 585 handpiece: 550 – 1100 nm ReMove 650 - 625 – 1100 nm
Pulse duration	sub-millisecond – 10ms
Pulse sequencing	2 – 6 pulses
Maximum Repetition rate	0.4 Hz
Maximum fluences	585 handpiece 40 J/cm <sup>2</sup> 650 handpiece 36 J/cm <sup>2</sup>
Spot sizes (ILG)	5x1 cm, 2.5x1 cm, 1x1 cm and 1 cm Diam Circular (optional)

---

#### System Specifications

#### ReMove LP

---

Wavelengths	1064 nm
Repetition rate	0.4 Hz
Pulse duration	Sub-millisecond - 5ms
Maximum fluence	300 J/cm <sup>2</sup>
Spot sizes	3, 6, 8 mm

---

#### System Specifications

#### ResurFACE®

---

Wavelengths	2940 nm
Repetition rate	0.4 Hz with Pulse Stacking
Treatment zone	10 mm diameter
Fractional output	40 microspots per cm <sup>2</sup>
Energy per microspot	Up to 20 mJ

---

#### System Specifications

#### ReMove Q

---

Wavelength	532 nm and 1064 nm
Fluence Range	0.6 - 17 J/cm <sup>2</sup>
Repetition rate	Up to 5 Hz
Pulse duration	12 ns
Spot sizes	3, 5 mm
Delivery	Q-switch handpiece

---

#### Physical Parameters

---

Dimensions	1235 x 763 x 520 mm (H x D x W)
Weight	84 kg (Approx.)
Mains supply	110 - 240 VAC, 12 A
Line frequency	50 - 60 Hz
Manufacturer	Lynton

\*Other voltages are available on request

Dynamix is a brand name of the device officially registered and certified as 3JUVE. The medical device CE-Certification under the European Medical Devices Directive (MDD) only relates to the medical indications for use detailed in the accompanying Operator Manual.

Call Lynton today on 01477 536 977  
or visit [www.lynton.co.uk](http://www.lynton.co.uk)



[Lynton Lasers Ltd](http://www.lynton.co.uk)

Lynton House, Manor Lane,  
Holmes Chapel, Cheshire, CW4 8AF

Tel: +44 (0)1477 536 977

[info@lynton.co.uk](mailto:info@lynton.co.uk)

[www.lynton.co.uk](http://www.lynton.co.uk)

CE  
0123

Aesthetics  
**AWARDS**  
Winner 2023